Simple Perfusion Kit

SPK, Simple Perfusion Kit

Low cost 4 or 8 Channel Perfusion Systems that come as a kit. The SPK-4 or SPK-8 are excellent choices for making simple perfusion systems. The kits are great for intracellular and extracellular recordings as well as patch clamp and imaging applications. They are especially useful for all types of bath perfusion applications and for the student lab.


SPK-4 & SPK-8 System Highlights

- 4 or 8 channel kits
- easy to assemble
- low cost
- easy adjustment of flow
- excellent for student setups
- · high quality parts used

The **SPK-4** or **SPK-8**, **Simple Perfusion Kit,** are low cost manual gravity fed perfusion systems. They are specifically designed for bath, tissue, cell, and Oocyte perfusion. They are excellent for use in electrophysiology and imaging studies as well as student labs. These kits include all of the parts necessary to make a four or eight channel perfusion system. The idea was to make a system that was flexible and low cost, yet professional. All parts come in separately marked bags






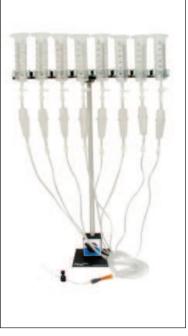

SPK-4 system

SPK-8 system

with clear instructions about how to assemble the kit. The kit includes a magnetic stand, a 4 / 8 channel reservoir bracket, 4 / 8 drip chambers, 4 / 8 manual 3 way luer valves, MMF-4 or MMF-8 Minimanifold™ output, PV-2 flow control valve, and misc. tubing and fittings. The flow valve allows the user to adjust flow at the output from 0ml/min up to 6ml/min. This valve equalizes the flow from any reservoir to the desired flow rate.

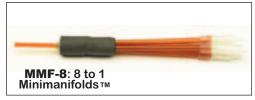

| SPK-4 and SPK-8 Ordering information: | | | |
|---------------------------------------|--|-------------|--|
| SPK-4 | 4 channel Simple Perfusion Kit - inc. magnetic stand, reservoir bracket, drip chambers, tubing, minimanifoldÔ, flow control valve, 3 way luer valves, accessories, instruction manual | SPK-8 | 8 channel Simple Perfusion Kit - inc. magnetic stand, reservoir bracket, drip chambers, tubing, minimanifoldÔ, flow control valve, 3 way luer valves, accessories, instruction manual |
| MMF-4 | 4 channel minimanifold™ outlet | MMF-8 | 8 channel minimanifoldÔ outlet |
| RESBRACK-4 | Replacement 4 channel reservoir bracket w/clamps for 60ml reservoirs | RESBRACK-8 | Replacement 8 channel reservoir bracket w/clamps for 60ml reservoirs |
| VM-STAND | Replacement magnetic stand w/2 poles | 60mISYRINGE | Box of 25 60ml reservoirs |
| SPK-DRIP | Replacement drip chambers - pack of 10 | 3WLV | 3 way stop cock luer valves |
| PV-2 | Flow control manual adjusted pinch valve | TUBING-2 | 1/16" ID x 1/8" OD tygon tubing - 10' |